

## Sterile SatPax® 670

Sterile Pre-wetted 55% Cellulose / 45% Polyester Nonwoven Cleanroom Wiper

Sterile SatPax® 670 combines Durx® 670 nonwoven wipers composed of a hydroentangled nonwoven blend of 55% cellulose and 45% polyester with a high saturation level of 70% IPA and 30% DI water.

This pre-wetted format provides a cost effective and easy to use solution versus traditional bulk handling of solvents, maintenance of squirt bottles and inconsistent wetting and cleaning associated with wetting a dry wiper.

Each package is then gamma irradiated sterile to a 10<sup>-6</sup> Sterility Assurance Level making the product ideal for cleaning critical sterile processing environments.



### Other Sterile Pre-wetted wipers

- Sterile SatPax® 1000
- Sterile SatPax®1200
- Sterile SatPax®-HA

### Key Attributes

- 55% cellulose / 45% polyester hydroentangled nonwoven blend
- No chemical binders in base material
- Pre-wetted with consistent 70% IPA/ 30% DI Water to a high saturation level
- Re-sealable solvent resistant packaging
- Wipers in c-folded configuration for single withdrawal
- Gamma irradiated and sterile validated to a 10<sup>-6</sup> Sterility Assurance Level per AAMI guidelines
- Certificate of Sterility available to download on the web by lot number **24/7**. Details: Expiration Date, Radiation Dosage, Lot Information, Meets cGMP requirements for traceability

### Benefits

- Low extractables and fiber and particle counts
- Smooth and durable with good wet strength
- Reduces alcohol usage and preparation / handling costs
- Reduces VOC emissions
- Increases cleaning efficiency
- Increases cleaning protocol consistency

### Environmental

- Reduces VOC (Volatile Organic Compound) emissions
- 55% cellulose fiber content is biodegradable
- Complies with RoHS
- Registered under REACH

### Applications

- Designed for use in ISO Class 5 and higher sterile cleanroom environments and USP <797> applications
- Designed for use in wet cleaning of critical surfaces where control of flammable solvents and flammable solvent concentrations is required
- Final cleaning of surfaces or products prior to manufacturing or packaging
- High saturation level is ideal for removing cleaning and disinfecting residues in sterile environments
- Meets AMS 3819C requirements, Class2, Grade A, Form 2

### Saturation Levels

The amount of solution contained in each wiper can be varied according to customer requirements. Higher saturation levels apply more solution to the surface during cleaning.

### Validated Sterile

ANSI/AAMI/ISO procedures have been followed for verification of irradiation dose with audits being performed to monitor bioburden levels and dose. The Sterility Assurance Level (SAL) or the probability of a viable organism being present on a product unit after sterilization is a SAL of 10<sup>-6</sup> (which is the chance of 1 non-sterile unit in 1,000,000).

[www.berkshire.eu.com](http://www.berkshire.eu.com)

Contact: Tel 44 (0) 870 757 2877  
 enquiries@berkshire.uk.com

Americas	Tel 1 413 528 2602	Fax 1 413 528 2614	info@berkshire.com
Europe	Tel 44 (0) 870 757 2877	Fax 44 (0) 870 757 2878	enquiries@berkshire.uk.com
SE Asia	Tel 65 6252 4313	Fax 65 6252 4312	enquiries@berkshire.com.sg
Japan	Tel 81 3 5827 2380	Fax 81 3 5827 2382	master@berkshire.co.jp

## Technical Data (In Dry State)

Attribute		Units	Value	Test Method
Basis Weight		g/m <sup>2</sup>	68.0	TAPPI T-410
Caliper		µm	264	TAPPI T-411
Fibers	≥100µm	fibers/cm <sup>2</sup>	160	IEST-RP.CC004.3, Sec 6.1.3 / Sec 6.2.2
Particles	≥0.5µm	x10 <sup>9</sup> /cm <sup>2</sup>	10	IEST-RP.CC004.3, Sec 6.1.3 / Sec 6.2.1
Sorbency	Capacity	mL/m <sup>2</sup>	320	IEST-RP.CC004.3, Sec 8.1 modified / Sec 8.2 modified
	Efficiency	mL/g	4.7	
	Rate	seconds	2	
Non-Volatile Residue	DI Water	g/m <sup>2</sup>	0.028	IEST-RP.CC004.3, Sec 7.1.2
	IPA	g/m <sup>2</sup>	0.0038	
Ions	Na <sup>+</sup>	ppm	62	IEST-RP.CC004.3, Sec 7.2.2
	K <sup>+</sup>	ppm	5.9	
	Ca <sup>++</sup>	ppm	22	
	Mg <sup>++</sup>	ppm	5.0	
	Cl <sup>-</sup>	ppm	31	

### Notes:

- Technical data represented in this table are typical values at the time of publication. These should not be used as product specifications.
- Due to differences in test methods applied and equipment utilized by different wiper manufacturers, valid product comparisons may only be obtained through side-by-side testing in the same test facility, under similar conditions
- Third party testing can be performed upon request

### Order Information:

Product	Number	Size	Shts/pk	Pks/cs	IPA/DI	Saturation	VOC % by Weight	Style
Sterile SatPax® 670	SSP670.001.24	9x9" (23x23cm)	30	24	70/30	60%	44%	C-fold

## Other Berkshire products




Wipers



Glove Liners



Mop Systems

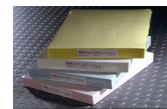

Documentation Systems



Face Masks



Swabs